DE Läusekamm Gebrauchsanweisung.....	2	IT Pettine anti pidocchi Istruzioni per l'uso	34
EN Lice comb Instructions for use.....	11	TR Bit Tarağı Kullanım Kilavuzu	42
FR Peigne anti-poux Mode d'emploi	18	RU Гребень от вшей Инструкция по эксплуатации.....	49
ES Lendreras Instrucciones de uso	26	PL Grzebień na wszy Instrukcja obsługi	57

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Read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

Table of contents

1. Included in delivery.....	11	7. Usage	15
2. Signs and symbols	11	8. Cleaning and maintenance.....	16
3. Intended use.....	12	9. Frequently asked questions.....	16
4. Warnings and safety notes	13	10. Disposal.....	16
5. Device description	15	11. Technical specifications	17
6. Initial use.....	15	12. Warranty/service.....	17

1. Included in delivery

Check that the exterior of the cardboard delivery packaging is intact and make sure that all contents are present. Before use, ensure that there is no visible damage to the device or accessories and that all packaging material has been removed. If you have any doubts, do not use the device and contact your retailer or the specified Customer Service address.

- 1 x HT 15 lice comb
- 1 x protective cap
- 1 x cleaning brush
- 1 x AA 1.5V battery
- 1 x set of these instructions for use

2. Signs and symbols

The following symbols are used in these instructions for use, on the packaging and on the type plate for the device and the accessories:

 WARNING	WARNING denotes danger with medium-level risk, which can lead to death or serious injury if ignored.
 CAUTION	CAUTION denotes danger with low-level risk, which can lead to minor or moderate injury if ignored.
	Product information Note on important information
	Observe the instructions Read the instructions before starting work and/or operating devices or machines
	Application part type B
	Do not dispose of batteries containing harmful substances with household waste

	Manufacturer
	Disposal in accordance with the Waste Electrical and Electronic Equipment EC Directive – WEEE
	Separate the packaging elements and dispose of them in accordance with local regulations.
	Marking to identify the packaging material. A = Material code, B = Material number: 1-7 = Plastics, 20-22 = Paper and cardboard
	Unique Device Identifier (UDI) for unique product identification
	Item number
	Batch designation
	Medical device
IP22	Protected against solid foreign objects 12.5 mm in diameter and larger, and against drops of water when the housing is angled up to 15°
Storage/Transport 	Permissible storage and transport temperature and humidity
Operating 	Permissible operating temperature and humidity
	CE labelling This product satisfies the requirements of the applicable European and national directives.
	Separate the product and packaging elements and dispose of them in accordance with local regulations.
	United Kingdom Conformity Assessed Mark
	Importer symbol

3. Intended use

Intended Purpose

The lice comb is designed to remove lice and/or nits when combing through the hair. When it comes into contact with a louse or nit, the lice comb emits a harmless electrical impulse that is imperceptible to humans, which kills the lice/nits. The lice comb is intended for domestic/private use only, not for commercial use.

Intended Users

The lice comb is only suitable for combing human hair of people of all ages.

Indications

The lice comb is used to kill and remove lice and/or nits in the event of a lice/nit infestation on a person's head.

Clinical benefit

The lice comb quickly and effectively kills and removes lice and nits without the use of chemical agents.

Contraindications

- Do not use the device on birthmarks or the scalp if the skin is broken.
- Do not use the device on persons with heart conditions, epilepsy or pacemakers.
- People with heart conditions, epilepsy or pacemakers must also not operate the device on other people.

Undesired side effects

- People with a sensitive scalp may experience slight tingling or itching during use. This is normal and not harmful to your health.
- Skin contact with the metal teeth during use may result in painful but harmless electrical impulses.
- The metal tines may cause injuries.

If used on several people:

- risk of transmitting live lice
- Cross infection of germs from a person being treated to another person

4. Warnings and safety notes

General warnings

- Use the device on dry, combed hair only.
- Do not use the comb if the metal tines are deformed as this may cause a short circuit in the comb.
- Keep packaging material away from children. Danger of suffocation!
- Do not allow children to play with the device.
- If the lice comb has been used on an animal, it must then never be used on humans.
- Store the device out of the reach of children.
- Children must always be supervised by an adult when using the device.
- After each use, clean the comb thoroughly using the cleaning brush – particularly before using the device on another person. For further information, see the chapter on cleaning and maintenance.
- If the lice comb is used by several people in succession without being cleaned properly, it may cause the transmission of diseases, such as mycoses (fungal disease).
- During use, ensure that the comb does not come into contact with the mouth, eyes or ears.
- Do not touch the metal teeth when the device is switched on
- This device can be used by children over the age of 12 and by people with reduced physical, sensory or mental skills or people lacking experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely and are aware of the consequent risks of use.
- Never operate the device during maintenance work.

- In order to prevent any damage to health, we would urgently advise against using the electric lice comb in the following cases without consulting a doctor:
 - With any implanted electrical devices (such as pacemakers)
 - If there are any metal implants. Make sure that no metallic objects such as necklaces can come in contact with the device during usage. If you wear any jewelry or piercings (e.g. ear piercing) near the application site, you must remove these before using the device.


General precautions

- Keep the device away from hot surfaces, and also protect it against direct sunlight.
- When not in use, always place the protective cap provided on the device.
- Operation in the vicinity of strong electromagnetic fields, for example next to a mobile phone, may cause malfunctions.
- Before use, inspect the device and accessories for visible damage. If you have any doubts, do not use the device and contact your retailer or the specified Customer Service address.
- Do not use the device near or in water.
- Never submerge the device in water or other liquids! Make sure that moisture does not get inside the device.
- Protect the device from impacts, moisture, dust, chemicals, marked temperature fluctuations and nearby sources of heat (ovens, heaters).
- Do not attempt to disassemble the device.
- Never place any heavy objects on the device.

Notes on handling batteries

- Always insert the batteries correctly, taking into account the polarity (+ / -). Keep batteries clean and dry and away from water. Always select the correct battery type.
- Never short-circuit batteries and battery compartment contacts.
- Never charge, forcibly discharge, heat, disassemble, deform, encapsulate or modify batteries.
- Never weld or solder on batteries.
- Never mix batteries of different manufacturers, capacities (new and used), size and type within a device.
- **Risk of explosion!** Failure to comply with the points mentioned above can result in personal injury, overheating, leakage, venting, breakage, explosion or fire.
- If a battery has leaked, put on protective gloves and clean the battery compartment with a dry cloth.
- If your skin or eyes come into contact with battery fluid, rinse the affected areas with water and seek medical assistance.
- **Choking hazard!** Keep batteries out of the reach of children. Seek medical attention immediately if swallowed.
- Never allow children to replace batteries without adult supervision.
- Store batteries away from metal objects in a well-ventilated, dry and cool place.
- Never expose batteries to direct sunlight or rain.
- Remove batteries from the device if it is not going to be used for a long period of time.
- Dispose of discharged batteries immediately and properly. Never dispose of batteries in fire.
- When disposing of batteries, keep batteries with different electrochemical systems separate.

5. Device description

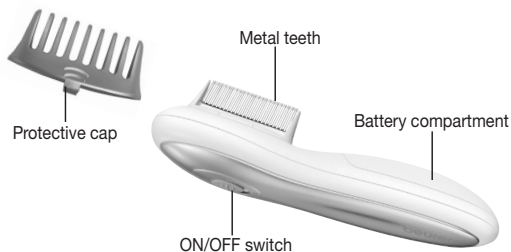6. Initial use

Inserting the battery

1. Slide the battery compartment lid open.
2. Insert the battery included in delivery into the battery compartment as shown inside the battery compartment. When inserting the battery, ensure that the correct polarity (+/-) is observed.
3. Then slide the battery compartment lid closed until you hear and feel it click into place.

7. Usage

Only use the lice comb on dry hair. Before use, comb the hair to be treated with a normal comb right to the tips until all tangles are removed.

1. Using your thumb and index finger, press on the notches on the side of the protective cap until it is released. Pull the protective cap up and off.
2. Move the ON/OFF switch to "1". A high-pitched beep sounds, which indicates that the device is switched on and ready for use.
3. Hold the lice comb at an angle of 45° to the scalp. Slowly move the lice comb from the roots to the tips of the hair. Slowly comb through all of the hair carefully and completely.
4. As soon as a louse comes into contact with the metal teeth of the lice comb, the beep audibly changes. If the high-pitched beep does not sound again, clean the metal teeth with the cleaning brush.
5. To end the application, move the ON/OFF switch to "0". The high-pitched beep is switched off.
6. Clean the lice comb after every use as described in the "Cleaning and maintenance" chapter.

NOTE

The lice comb removes hatched lice as well as unhatched nits (the eggs of the lice). Repeat the application every day for two weeks so that the lice are completely eliminated. In order to prevent lice returning, we recommend repeating the application once or twice a week for a further month. This means that any remaining lice can be removed.

8. Cleaning and maintenance

CAUTION

- Remove the protective cap for better cleaning of the metal teeth. Carefully clean the metal teeth after every use using the cleaning brush provided. Check the lice comb for any remnants between the metal teeth.
- Clean the device (including the metal teeth) after use with a soft, slightly damp cloth. If it is very dirty, you can also moisten the cloth with a mild soapy solution. Do not use any chemical or abrasive cleaning products.
- Ensure that no water enters the device. If this should happen, remove the batteries and only use the device again when it is completely dry.

9. Frequently asked questions

Who can be affected by head lice?

Anyone can be affected by head lice, not only children. Having a case of lice is also not caused by a lack of hygiene.

How do you get head lice?

Lice cannot fly or jump; they scuttle and crawl. Lice are transferred by direct contact with a person that already has lice (for example when two heads come into contact).

Where do head lice come from?

School is usually not the source of transfer. Lice are found much more often in communities outside of school.

Who should be checked and how?

If you have been informed of lice in your community, check all members of your family at least once a week and inform your friends.

Why is it important to ascertain whether lice are present immediately?

If you find out promptly that lice are present, it is easier to remove them as their numbers will be lower. There is also a lower probability that the lice will be transferred to other family members or acquaintances.

10. Disposal

For environmental reasons, do not dispose of the device in the household waste at the end of its useful life. Dispose of the device at a suitable local collection or recycling point in your country. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment).

If you have any questions, please contact the local authorities responsible for waste disposal. Batteries must not be disposed of in the household waste. They may contain poisonous heavy metals and are subject to special refuse treatment. The codes below are printed on batteries containing harmful substances:

Pb = battery contains lead,

Cd = battery contains cadmium,

Hg = battery contains mercury.

11. Technical specifications

Model:	HT 15
Battery:	1 x 1.5V AA batteries
Weight:	Approx. 52 g (without batteries)
Dimensions:	158 x 37 x 50 mm
Permissible operating temperature and humidity:	10°C to 40°C, 30% to 85% relative humidity
Expected service life of the device	Information on the service life of the product can be found on the homepage

The batch number is located on the device or in the battery compartment.
Subject to technical changes.

The device complies with well as the European standard EN 60601-1-2 (Group 1, Class B in accordance with IEC 61000-4-2, IEC 61000-4-3) and is subject to particular precautions with regard to electromagnetic compatibility. The device is suitable for use in all environments listed in these instructions for use, including domestic environments.

Notes on electromagnetic compatibility

- The device is suitable for use in all environments listed in these instructions for use, including domestic environments.
- The use of the device may be limited in the presence of electromagnetic disturbances. This could result in issues such as error messages or the failure of the display/device.
- Avoid using this device directly next to other devices or stacked on top of other devices, as this could lead to faulty operation. If, however, it is necessary to use the device in the manner stated, this device as well as the other devices must be monitored to ensure they are working properly.
- The use of accessories other than those specified or provided by the manufacturer of this device can lead to an increase in electromagnetic emissions or a decrease in the device's electromagnetic immunity; this can result in faulty operation.
- Keep portable RF communication devices (including peripheral equipment, such as antenna cables or external antennas) at least 30 cm away from all device parts, including all cables included in delivery.
- Failure to comply with the above can impair the performance of the device.

12. Warranty/service

Further information on the warranty and warranty conditions can be found in the warranty leaflet supplied.

Notification of incidents

For users/patients in the European Union and identical regulation systems (EU Medical Device Regulation (MDR) 2017/745), the following applies: If during or through use of the product a major incident occurs, notify the manufacturer and/or their representative of this as well as the respective national authority of the member state in which the user/patient is located.

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UK-Importer: Beurer UK Ltd.
Beurer GmbH • Söflinger Str. 218 • 89077 Ulm www.beurer.com