A E N O

Quick Start Guide


Smart Toothbrush DB1S / DB2S


(Models: ADB0001S,
ADB0002S)

27.08.2024 | v 2.2.0


ENG.....	8
ARA.....	13
BOS.....	16
BUL.....	20
CES.....	25
DEU.....	29
ELL.....	34
EST.....	38
FRA.....	42
HRV.....	47
HUN.....	51
HYE.....	56
ITA.....	60
KAT.....	65
KAZ.....	70
LAV.....	75
LIT.....	79
NLD.....	84
POL.....	89
POR.....	94
RON.....	99
RUS.....	104
SLK.....	109
SLV.....	113
SPA.....	117
SRP.....	122
UKR.....	126
UZB.....	131

A

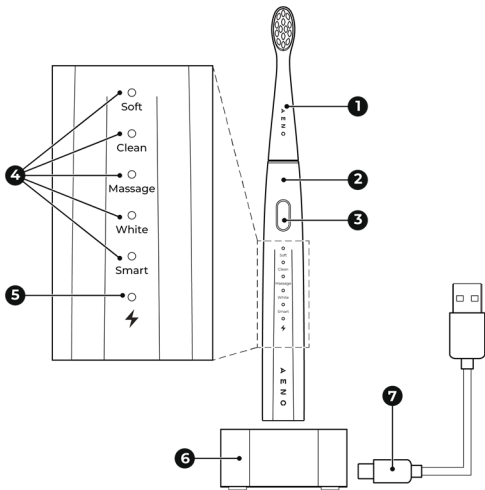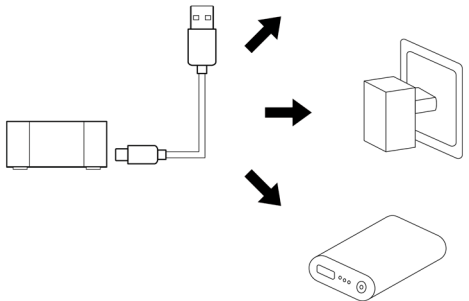
B

C

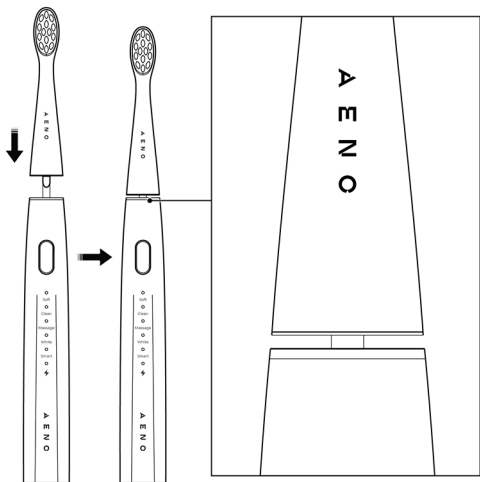D

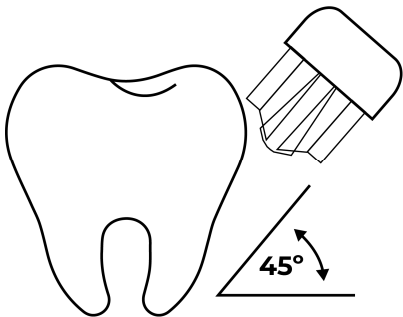ENG AENO™ Smart Toothbrush is intended for daily oral hygiene.

Models: ADB0001S (White), ADB0002S (Black).

Technical Specifications

Toothbrush. Input: 5.0 V (DC), 0.3 A, 1.5 W (max.). Battery: non-removable, lithium-ion, 1,400 mAh, 3.7 V, full charge time – up to 5.5 h. Operating time after a full charge: up to 90 days. Vibration intensity: up to 46,000 vibrations per minute. Ingress protection rating: IPX7. Noise level (at a distance of 0.3 m): ≤ 60 dB. Number of operating modes: 12*. Control: manual, remote. BT 4.2 (2402–2480 MHz). Transmission power: +4 dBm / 2.5 mW. Dimensions of the toothbrush with the brush head (diameter × height): 25.8×240 mm. Casing material: ABS plastic. Stiffness of the bristles: medium. Service life of brush heads: 3 months. Weight: 105 g. Operating and storage conditions: temp. 0...+40 °C, RH 45–95 %.

Charging Station. Input: 5.0 V (DC), 1.0 A, 5.0 W (max.). Output: 5.0 V (DC), 0.3 A, 1.5 W (max.). Ingress protection rating: IPX5. Material: ABS plastic. Dimensions (diameter × height): 45×21 mm. Weight: 20 g.

Package Contents (see Fig. A)

AENO™ Smart Toothbrush, replacement brush heads (3 pcs.), charging station, USB (A) / USB (C) cable (0.8 m), Quick Start Guide.

Note. Power adapter is not included. The output parameters of the adapter must be 5.0 V (DC), 1.0 A.

Device Elements and Accessories (see Fig. A)

1 – replacement brush head, 2 – handle, 3 – button, 4 – operating mode indicators, 5 – battery charge indicator, 6 – charging station, 7 – USB (A) / USB (C) cable.

If you experience any questions or difficulties in using the AENO™ device, please email support at **support@aeno.com** or chat with us online at **aeno.com/service-and-warranty**. Our specialists will help you solve the problem, and you will not have to waste time and effort visiting the store.

* 4 operating modes are available on the device itself, and 8 operating modes are available in the AENO app.

Restrictions and Warnings

ATTENTION! If you have a pacemaker, oral diseases, or if you have just undergone oral surgery, consult your doctor before using the device.

ATTENTION! To extend battery life, do not store the device on a charging station connected to a power source.

This device is not intended for use by children under the age of 6 or by persons with reduced physical, mental or intellectual abilities unless they have sufficient experience and knowledge of how to use the device and are supervised by a person responsible for their safety. Do not allow children to play with the device and its accessories. Do not drop or throw the device. Do not turn on or use a damaged device. If your gums begin to bleed after using the device and the bleeding continues for several days, contact your dentist. Use only the original charging station to charge the toothbrush. Make sure that the voltage and frequency specified in the technical documentation correspond to the power supply parameters. Each family member should use an individual brush head. Replace the brush head at least once every 3 months. Use only original replacement brush heads. If you are not using the device for a long time, charge it every 2-3 months to prevent deep discharge of the battery.

Statuses of the Battery Charge Indicator

Color	Status	Meaning
Green	On	The device is turned on, the battery level is medium/high
		The device is placed on a charging station that is connected to a power source, the battery is fully charged
	Flickering	The battery is charging
Red	On	The device is turned on, the battery level is low
		The battery is fully discharged, the device does not turn on
-	Off	The device is turned on or placed on a charging station that is not connected to a power source

States of the Operating Mode Indicators

You can tell the current operating mode of the toothbrush by the solid white indicator located above the name of the mode.

Button Functions

Action	Meaning
Press once (device is turned off)	Turning the device on
Press once (it has been less than 5 seconds since the device was turned on)	Switch operating modes
Press once (it has been more than 5 seconds since the device was turned on)	Turn the device off

Operating Modes

- “Soft” for delicate, gentle brushing of the teeth;
- “Clean” for daily brushing of the teeth;
- “Massage” for massaging the gums;
- “White” for teeth whitening;
- “Smart” for operating the device according to the scenario defined in the AENO app.

Preparing for Use

Charging. Carefully open the package and take out the device and its accessories. Clean them (see “Cleaning and Maintenance”). Place the charging station (6) on a dry, horizontal surface and connect it to a power source (see Fig. B). Place the handle (2) on the charging station. The battery charge indicator (5) will start to flicker green. When the indicator turns solid green, the battery is fully charged and the device is ready for use. When charging is complete, disconnect the charging station from the power source.

Installing the Brush Head. Remove the handle from the charging station and install the brush head (1) as shown in Figure C.

Choosing the Operating Mode. Turn on the device by pressing the button (3). Continue to press the button to switch operating modes. The operating mode

can be changed within 5 seconds after the device was turned on. After selecting an operating mode, press the button again to turn off the device.

Note. The last selected operating mode is saved in the device memory and will be activated the next time the device is turned on.

Using the Device

Wet the bristles of the brush head and apply the desired amount of a dental care product to the bristles. Point the bristles of the brush head toward the tooth surface at an angle of about 45° (see Fig. **D**). Turn on the device by pressing the button once. Brush your teeth without putting too much pressure on the tooth enamel. Every 30 seconds, the device will notify you to change the brushing zone by briefly stopping the vibration. Change the brushing zone after each notification. After 2 minutes of operation, the device will turn off automatically.

Cleaning and Maintenance

ATTENTION! Turn off the device and disconnect the charging station from the power source before cleaning the device and its accessories.

ATTENTION! Do not wash the device or its accessories in a dishwasher.

ATTENTION! Do not use chemical and aggressive detergents, abrasive pastes, products containing acids and solvents or metal sponges to clean the device and its accessories.

Cleaning the Brush Head. Rinse the brush head, including the bristles, under running water.

Cleaning the Handle. Rinse the handle under running water or wipe it with a soft cloth dampened with water. Then wipe dry.

Cleaning the Charging Station. Wipe the station with a soft cloth dampened with water. Then wipe dry.

Connecting to AENO Mobile App

- Download the AENO mobile app from Google Play or App Store to your smartphone.
- Register a new account or log in to an existing account.

- Press the button on the device to turn it on. Make sure that the battery of the device is charged. The battery charge indicator should be green. Press the power button again to turn off the device.
- Place your smartphone as close to the device as possible.
- Tap “+” in the top right corner of the app's home screen, then “Add Device”.
Note. If you receive system notifications to turn on any features on your smartphone, follow the recommendations in the notifications.
- The device will be found automatically and its image will appear on the screen. Tap the “Add” button.
- In the next window, tap “+” next to the toothbrush image to add it to your account. Wait for the message “Successfully added” to appear and tap “Done”.

A detailed description of the device, its operating modes, as well as information on how to use the device in the app can be found in the complete User Manual available at aeno.com/documents.

Troubleshooting

The toothbrush does not turn on. Possible cause: the battery is discharged. Solution: charge the device (charging time – up to 5.5 hours).

The toothbrush is not charging. Possible causes: the charging station is not connected to a power source; the USB cable is damaged. Solution: connect the charging station to a power source (see Fig. **B**); check the integrity of the USB cable.

The app cannot find the toothbrush. Possible cause: BT is turned off on the smartphone. Solution: make sure BT is activated on the smartphone.

ATTENTION! If none of the suggested solutions helped you eliminate the problem, contact the supplier or an authorized service center. Do not disassemble the device or attempt to repair it yourself.

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ENG ASBISC reserves the right to modify the device and to make changes and additions to this document without prior notice to users.

The warranty period and service life is 2 years from the date of retail sale of the device.

Manufacturer: ASBISC Enterprises PLC, Iapetou, 1, Agios Athanasios, 4101, Limassol, Cyprus. Made in China. All trademarks and brands mentioned herein are the property of their respective owners. The date of manufacture is indicated on the package.

Up-to-date information and a detailed description of the device, as well as connection instructions, certificates, information on companies that accept quality claims and warranties is available at aeno.com/documents.

SIMPLIFIED DECLARATION OF CONFORMITY. Hereby, ASBISC declares that the device is in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following internet address: aeno.com/userfiles/files/Certificates/CE/DB1S-DB2S.pdf.

RoHS Directive. The device complies with the requirements of the Directive 2011/65/EU on the Restriction of Hazardous Substances (RoHS), including the requirements of the RoHS Directive 2015/863/EU.

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ENG Recycling information. These symbols indicate that you must follow the Waste Electrical and Electronic Equipment (WEEE) and Waste Battery and Battery Regulations when disposing of the device, its batteries and accumulators, and its electrical and electronic accessories. According to the rules, this equipment must be disposed of separately at the end of its service life. Do not dispose of the device, its batteries and accumulators, or its electrical and electronic accessories together with unsorted municipal waste, as this will harm the environment. To dispose of this equipment, it must be returned to the point of sale or turned in to a local recycling facility. You should contact your local household waste disposal service for details.

For devices with a non-removable battery: under normal or abnormal operating conditions, as well as in the event of a single malfunction, the charging voltage and current must not exceed their maximum permissible values. If the battery temperature exceeds the maximum allowable temperature, then, under abnormal conditions, the device will stop charging. If the battery temperature drops below the minimum allowable temperature during charging, then, under abnormal conditions, the current will be limited.

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